

## Allies











## Connect

- alien\_innovations
- in ALIEN INNOVATIONS
- alieninnovations.org
  - connect@alieninnovations.org



USB-C port at the back for charging



Turn off the device while charging



Use the cables provided with the device



10-15 mins of full charge time




Don't use adapters of more than 12amp



## MAVINAUT

Empowering visually impaired enhancing mobility

- NAVINAUT is a torch like device for visually impaired for ease of their mobility without a stick
- It detects obstacles while walking and vibrates to alert them
- The torch type compact design makes it easy to carry
- The device is ergonomically designed for ease of use
- The compact design also helps in detecting smaller objects
- As it is a hand device it helps them detect obstacles in all directions
- This devices ensures peoples safety without the need of carrying a stick

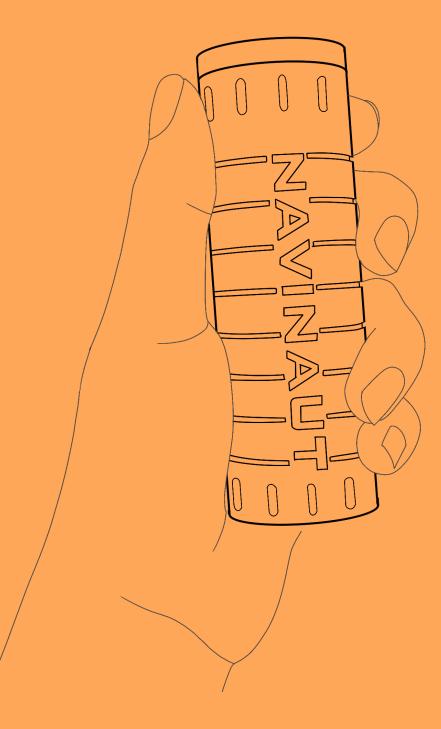



To turn on the device, press the button at the back side of the device



Once the button is pressed the device starts detecting the obstacles



To turn of the device press the same button